
AeroChamber^{Plus} Z STATTM

Anti-Static Valved Holding Chamber

MOUTHPIECE
LARGE MASK

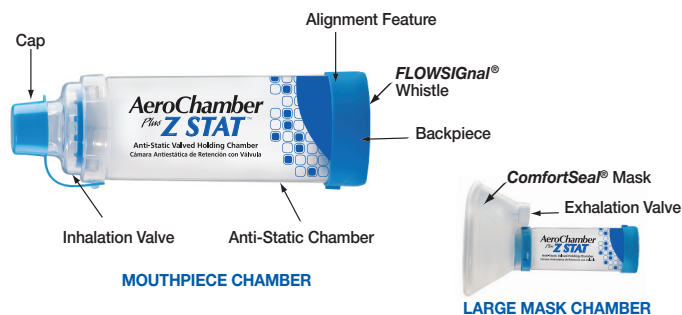INDICATIONS FOR USE

This product is intended to be used by patients who are under the care or treatment of a physician or licensed healthcare professional. The device is intended to be used by these patients to administer aerosolized medication from most pressurized Metered Dose Inhalers. The intended environments for use include the home, hospitals and clinics.

INSTRUCTIONS FOR USE

THIS PRODUCT CAN BE USED RIGHT OUT OF PACKAGE. BEFORE USE, ENSURE THESE INSTRUCTIONS AND THE INSTRUCTIONS SUPPLIED WITH THE METERED DOSE INHALER (MDI) HAVE BEEN READ.

MOUTHPIECE	1	2	3	4	5	6	7
LARGE MASK	1	2	3	4	5	6	7

CLEANING INSTRUCTIONS FOR MASK AND MOUTHPIECE CHAMBERS

THIS PRODUCT CAN BE USED RIGHT OUT OF THE PACKAGE AND THEN CLEANED WEEKLY.

1	2	3	4

Notes:

- Storage and operating range 5° C - 40° C (41° F - 104° F) at 15 to 95% relative humidity.
- Product may need to be replaced after 12 months of use. Environmental conditions, storage and proper cleaning can affect product life span.
- Do not share this medical device.
- Clarity of the chamber is a result of the properties of the StatBan® anti-static material.
- If you notice medication build-up in your chamber, wash the inside of the chamber gently with a soft cloth.
- Dishwashing with overly dirty dishes is not recommended.
- Dishwasher validated up to 158°F (70°C).

CAUTIONS:

- PRODUCT MAY BE PERMANENTLY DAMAGED IF BOILED, STERILIZED OR CLEANED IN A DISHWASHER AT TEMPERATURE ABOVE 158°F (70°C).
- Do not leave the chamber unattended with children.
- Federal (USA) law restricts the sale of this device on or by the order of a physician.

Product Reorder Number:
MOUTHPIECE: 79710/case 10, 79750/case 50
LARGE MASK: 80710/case 10

Rx Only

This device has not been made (or manufactured) with BPA (Bisphenol A) or natural rubber latex.

Manufactured by:

monaghan®

Monaghan Medical Corporation
Plattsburgh, NY 12901
www.monaghanmed.com

© 2013 Monaghan Medical Corporation™ and ® are trademarks and registered trademarks of Monaghan Medical Corporation or an affiliate of Monaghan Medical Corporation. This device was Made in the USA from US and Canadian parts.

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