



## STRIVE® DUAL ZONE PEAK FLOW METER



STRIVE® DUAL ZONE PEAK FLOW METER (PFM) provides accurate measurements for patient populations from 5 years old through adult. With just one device, the STRIVE® DUAL ZONE PFM can accommodate measurements from 50-800 L/min. Designed for individual or multi-patient use, the STRIVE® DUAL ZONE PFM is easy to clean for high-level disinfection and has optional disposable mouthpieces with one-way valves, for clinical applications.

### The Monaghan Advantage

- ❖ Innovative, cost effective, small and portable device makes it easy to track day-to-day changes in breathing patterns.
- ❖ Provides accurate measurements for patient populations from 5 years old through adult.
- ❖ Measurements from 50-800 L/min accommodated.
- ❖ Suitable for multi-patient use with optional disposable mouthpieces with one-way valves.
- ❖ Meets the most recent technical Standards of Spirometry for portable peak flow meters established by the American Thoracic Society.

MEETS U.S. AND  
INTERNATIONAL  
TECHNICAL  
STANDARDS

### THE FACTS

**The Strive® PFM - the one device that can simplify inventory needs and reduce costs.**

**Meets the most recent technical Standards of Spirometry for portable peak flow meters established by the American Thoracic Society.<sup>1</sup>**

**Meets International Standard for peak expiratory intended for the proper assessment of pulmonary function.<sup>2</sup>**

#### REFERENCES

1. ATS/ERS Task Force: "Standardization of Lung Function Testing," 2005.
2. International Standard ISO 23747:2015: "Anaesthetic and Respiratory Equipment - Peak Expiratory Flow Metrics for Assessment of Pulmonary Function in Spontaneously Breathing Humans."

518.561.7330 phone  
800.833.9653 toll free  
518.561.5088 fax

Monaghan Medical Corporation  
Plattsburgh, New York 12901  
[monaghanmed.com](http://monaghanmed.com)



## THE FACTS (continued)


**The Strive® Dual Zone Peak Flow Meter meets the educational guidelines recommended by the National Asthma Educational Prevention Program.**

**The Strive® PFM combines a low flow scale and standard scale in one device, providing the option to replace two devices with one device, reducing the number of items to inventory. With only one device to choose from and one device to inventory, the Strive® PFM saves both time and money.**

Accommodates measurements in two distinct zones: 50-350 L/min and 200-800 L/min

Incorporates the NAEPP three color zones and is equipped with a low-flow and high-flow scale

Adjustable, color coded zone markers can be moved to delineate the patient's zones.

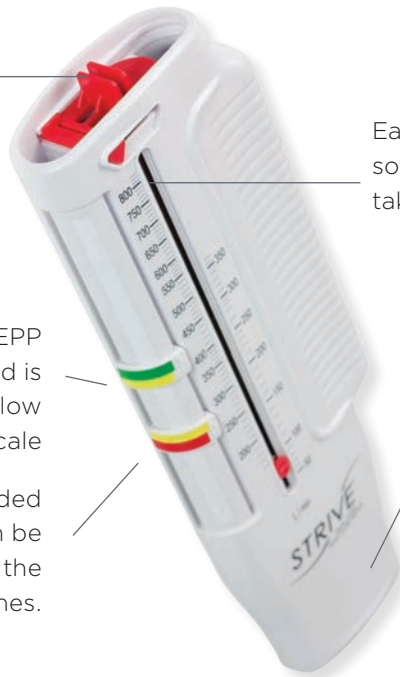

Easy to see numbers so readings can be taken quickly

Ergonomically shaped mouth-piece assists in creating a secure seal.

Multi-patient use with optional disposable mouthpieces.

## CLEANING IS EASY

- Designed for one year of use in the home, clinic or hospital environment
- Approved for the top shelf dishwasher cycle and high level disinfectants, such as Cidex<sup>+</sup> or Revital-OX<sup>+</sup> Resert<sup>+</sup>.

## ORDERING INFORMATION

PART #	NAME/DESCRIPTION	UNITS PER CASE
96610	STRIVE® Dual Zone Peak Flow Meter (PFM)	cs/10
90200	Disposable mouthpiece	cs/200

Contact your Monaghan Sales Territory Manager for more information or a demonstration.

